


EKG 12 lead placement

Electrode	Placement
V1	4 th Intercostal space to the right of the sternum
V2	4 th Intercostal space to the left of the sternum
V3	Midway between V2 and V4
V4	5 th intercostal space at the midclavicular line
V5	Anterior axillary line at the same level as V4
V6	Midaxillary line at the same level as V4 and V5
RL	Anywhere above the right ankle and below the torso
RA	Anywhere between the right shoulder and the wrist
LL	Anywhere above the left ankle and below the torso
LA	Anywhere between the left shoulder and the wrist

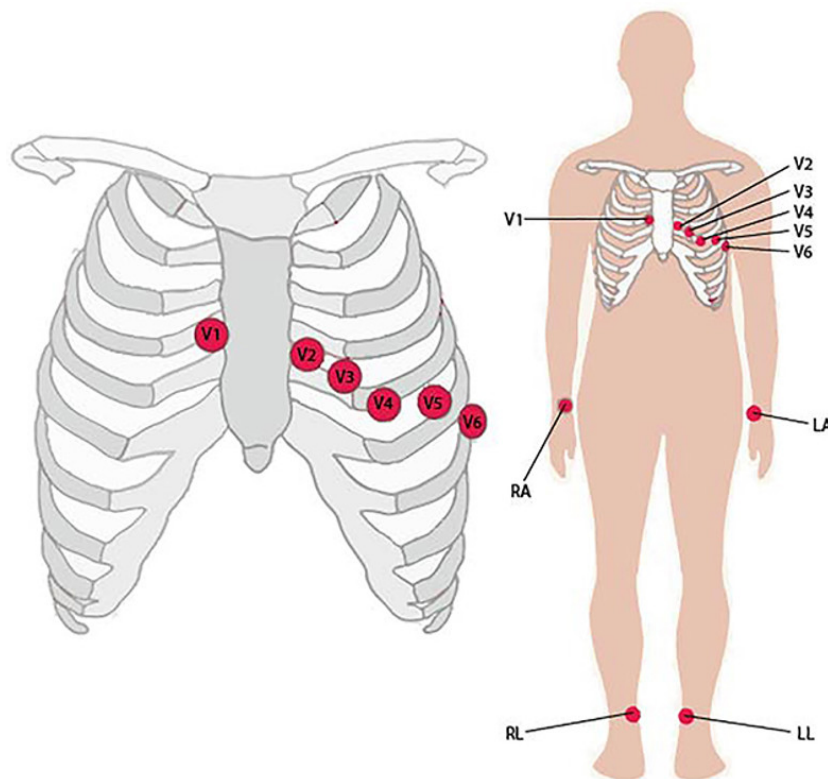

Created 5/2024